



MRC PPU  
**Reagents  
and Services**

## **ANTIBODY TESTING RESULTS** **Standard Reporting Template**

**Antibody Name:** Anti-IAV-PA

**Full Antigen Name:** PA (polymerase A) protein of influenza A virus (IAV) strain A/Puerto Rico/8/1934(H1N1) ('PR8')

**Antigen Species:** Virus

**Antigen:** GST-IAV-PA (1 – 716) DU70818, recombinantly expressed in Bacteria

**Purification:** serum heated treated, filtered and affinity purified against MBP-IAV-PA (1 - 716) DU70826

**Sheep Number:** DA186

**Bleeds Tested:** 1<sup>st</sup>, 2<sup>nd</sup>, 3<sup>rd</sup>, 4<sup>th</sup>, 5<sup>th</sup>

**Recommended Bleed:** All bleeds suitable for applications listed.

## Immunoblotting:

### Method

MDCK cells were infected with IAV strains PR8 or A/WSN/33 ('WSN') at an MOI of 3 PFU/cell, or mock infected.

293T cells were transfected with PR8 pDUAL:PA (a bidirectional reverse genetics plasmid that encodes the PA protein) or mock transfected, using Lipofectamine 2000 (Thermo Fisher) using cell density and plasmid mass as recommended in the manufacturer's protocol.

Cells were lysed in Laemmli buffer at 24 h post-infection or 48 h post-transfection, separated by electrophoresis on 12 % polyacrylamide gels and transferred to nitrocellulose membranes.


Membranes were blocked with 5 % milk in PBS-T (phosphate buffered saline with 0.1 % Tween 20) overnight at 4 °C, rinsed in PBS-T then probed with primary antibodies at 1 µg/ml in PBS-T for 1 h at room temperature. Membranes were washed ×3 in PBS-T followed by further incubation with anti-sheep IR680 or IR800 (Thermo Fisher) at 1 in 10,000 in PBS-T, for 1 h at room temperature.

Membranes were scanned with a Licor Odyssey CLx Infrared imaging system.

### Results


IAV PA = 82.60kDa, PA-N155 = 64.49kDa, PA-N182 = 61.27kDa, PA-x = 29.43kDa.. Ladder is the Page Ruler Prestained NIR Protein Ladder (Thermo Fisher).

#### (1) Infected cells

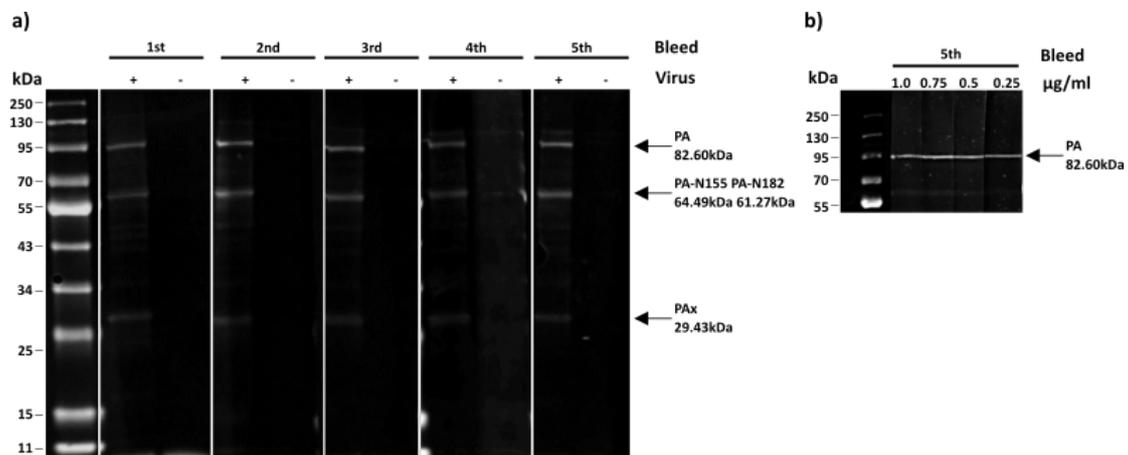

- Strong and specific signal was detected for all bleeds at 1 µg/ml. Antisera also bind to PA-N155/PA-N182 and PA-X proteins.
- Dilution series using 5<sup>th</sup> bleed; 0.25 µg/ml up to 1 µg/ml.

#### (2) Transfected cells.

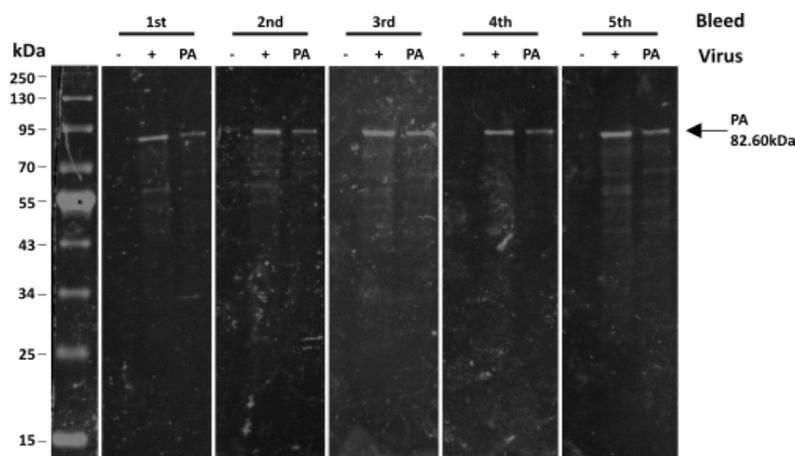

Strong and specific signal detected using all bleeds at 1 µg/ml.

**Recommendation**

All antisera can be used at 0.25  $\mu\text{g/ml}$  to 1  $\mu\text{g/ml}$  for western blot with infected cell lysates and at 1  $\mu\text{g/ml}$  for western blot with transfected cells.

## Immunofluorescence:

### Method

A549 and 293T cells were seeded at  $2 \times 10^5$  cells/well in 24 well plates containing 13 mm coverslips.

For 293T cells, coverslips were first coated with poly-D-Lysine to promote adhesion.

The day after seeding, A549 cells were infected with PR8 at an MOI of 2 PFU/cell and 293T cells were transfected with as above using Lipofectamine 2000, following the manufacturer's protocol.

Cells were fixed at 24 h post-infection or post-transfection in 4 % formaldehyde in PBS, then permeabilised using 0.2 % Triton-X100 in PBS/2 % fetal bovine serum (FBS) before blocking for up to 1 h in PBS/2 % FBS.

All antibody dilutions were carried out in PBS/2 % FBS. Cells were incubated with primary antibodies at 1:50 for 1 h at room temperature, washed 3 times in PBS/2 % FBS then probed with donkey anti-sheep IgG (H+L) Cross-Adsorbed Secondary antibody, AlexaFluor 647 (Invitrogen) at 1:500 and DAPI at 1:1000 for 45 minutes at room temperature.


Cells were washed 3 times in PBS/2 % FBS and once in PBS, excess moisture was removed before mounting coverslips on slides using ProLong Gold Antifade mounting agent (Thermo Fisher).

Mounting agent was allowed to harden overnight before observation using a Zeiss 710 confocal microscope.


### Results

A strong, specific signal was detected for all bleeds in both infection and transfection. Images show 5<sup>th</sup> bleed antibody at 24 h post-infection and 4<sup>th</sup> bleed antibody at 48 h post-transfection; similar results were observed for the other bleeds.

a)

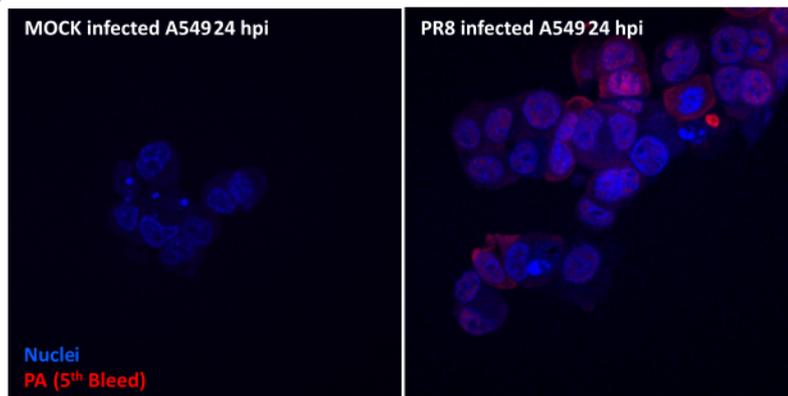

b)

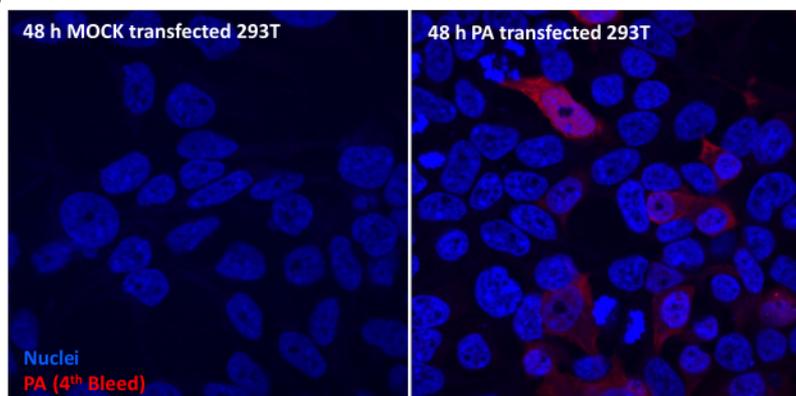

- a) Strong signal detected using all bleeds at 1 in 50 dilution in PR8 infected cells.
- b) Strong signal detected using all bleeds at 1 in 50 dilution in PA transfected cells.

### Recommendation

All 5 bleeds are suitable for use in IF of both infected and transfected cells at a 1 in 50 dilution.

## ELISA:

### Method

96 well ELISA plates were coated overnight at room temperature with lysates of PR8 infected or uninfected MDCK cells, with loading amounts normalised by tubulin concentration.

Plates were blocked for 2 h with 5 % milk in PBS-T.

All antibody dilutions were carried out using PBS-T.

Dilution series of the test antisera were added to the plates, prepared in duplicate from an initial concentration of 50 µg/ml using three-fold serial dilutions. Plates were incubated for 2 h at room temperature and washed 3 times in PBS-T.

Donkey Anti-Sheep IgG (H+L) Horseradish Peroxidase (HRP) Conjugate (Life Technologies/Novex) was added to wells at 1 in 20,000 dilution and incubated for 1 h at room temperature.


OPD substrate (Merck) was then added to wells and left to develop for 15 minutes in the dark and at room temperature, before stopping the reaction with 10 % H<sub>2</sub>SO<sub>4</sub>.

Absorbance (A490) was measured immediately using a Varioskan Lux (Thermo).

### Results:

Results were plotted as average absorbance at 490nm at each antibody concentration with results obtained using uninfected lysate coated plates subtracted from results obtained using infected lysate coated plates to correct for any non-specific binding.

Little/no signal above background was detected using any of 1<sup>st</sup> to 5<sup>th</sup> bleeds anti-PR8-PA using this method.

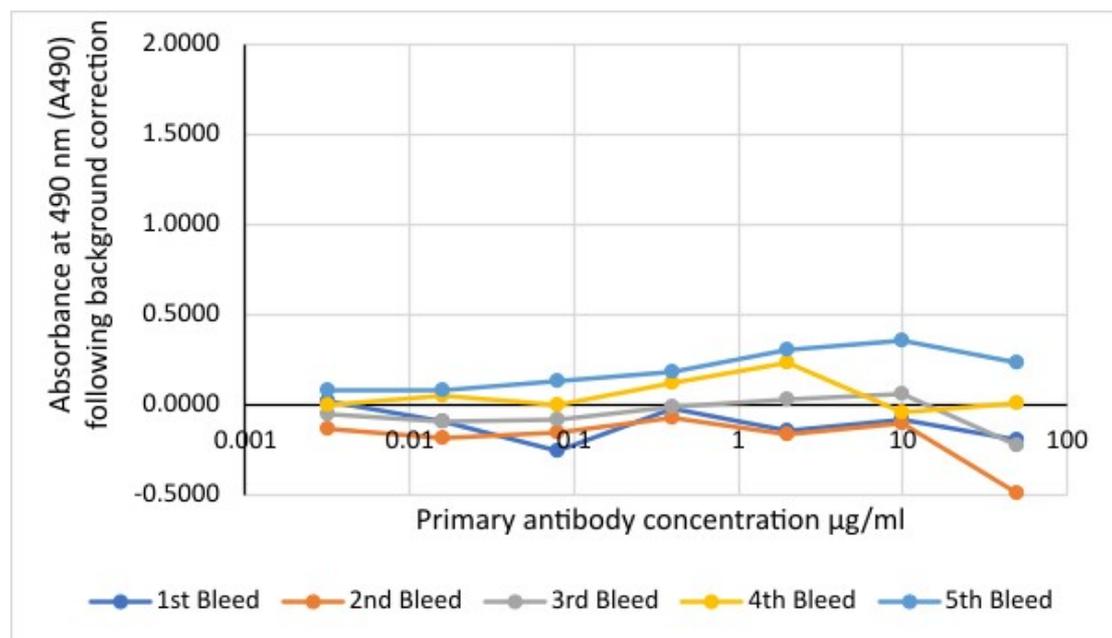

### Recommendation

Antisera not suitable for use in ELISA using this method

### PUBLICATIONS:

None to date

### Generated by:

Name: Sarah Cole      Date: 04/09/2023      Group Leader: Edward Hutchinson